



HETTINFO II

Documentation system
for blood banks



HETTINFO II DOCUMENTATION SYSTEM

— Transparency and traceability in blood banks

HettInfo II documents important work steps before and during centrifugation – this ensures transparency and traceable processes. The documentation system is operated via a touch display, which guides the user step by step through the centrifugation and records all significant process data.

All documented data points are temporarily stored locally and can be exported directly after the centrifugation or at a defined time as a CSV file. The data can be stored on the network or exported to a USB stick. This makes HettInfo II system independent and facilitates further processing into another software system.




— Good reasons to choose HettInfo II:

- 1 User-friendly touch display**
HettInfo II can be operated simply and easily via the touch display – even with laboratory gloves.
- 2 Reduced error rate**
With the step-by-step instructions, the user always knows exactly what to do next. In the event of incorrect operation or errors, a visual and audible message is issued immediately.
- 3 Automatic data export**
Data obtained can be exported after each run or in specific times.
- 4 Easy integration**
The exported CSV file is easy to process and work independently for greater compatibility with an existing system.



DOCUMENTATION PROCESS

- The HettInfo II process at a glance

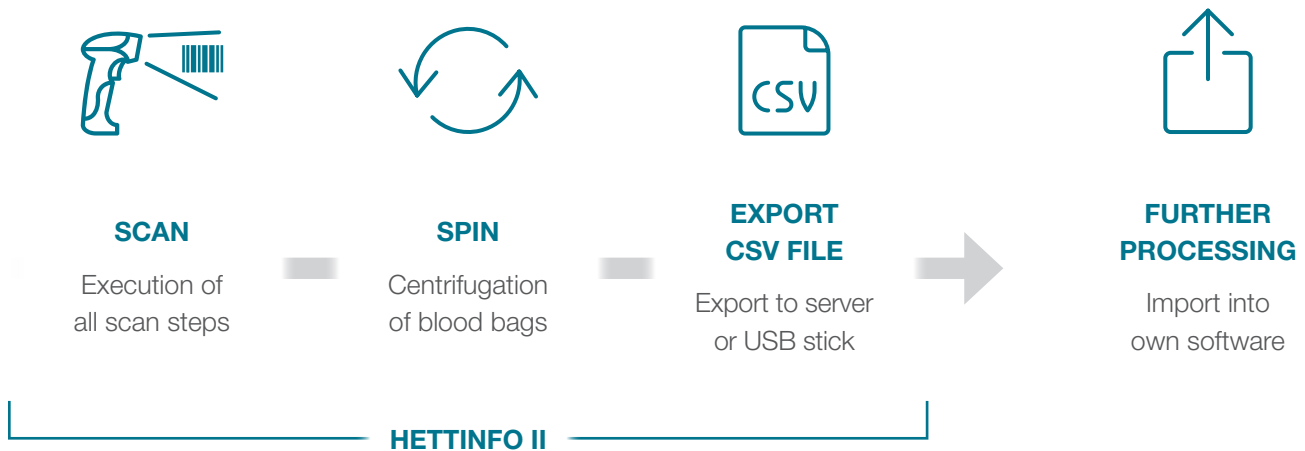

TECHNICAL FEATURES

- All technical features of HettInfo II

4.3" touch display	Process status display, error messages, settings
User	Any number of users possible
Data logging Scan Module	User ID, blood bag ID, program number
Run data logging	Centrifuge names, centrifuge data, date, start and end time, temperature, radius, acceleration time, braking time, total running time, RPM, RCF, integral RCF, cycle time, brake deactivation speed, end result
File format	CSV
Ports	LAN, USB
Backup	Redundant data backup
Logbook	Documentation of error messages. Export as text file on USB stick possible.



MODULES

— HettInfo II

With HettInfo II you can record running data of your centrifuge. To record a complete centrifugation process, we recommend HettInfo II plus Scan Module.



Articles	Includes
HettInfo II	Firmware, touch display

Cat. No.
0955

— Scan Module

The bar-code scanner included in the Scan Module records all data reliably. The ergonomic holder ensures easy handling and offers secure protection when not in use.

Bar-code scanners from the HettInfo I series are compatible with HettInfo II and can still be used.



Articles	Includes
Scan Module	Bar-code scanner, bracket, cable clamps

Cat. No.
0959

— Add-on Kit

With the Add-on Kit, ROTO SILENTA 630 RS models from model year 01/2008* can be upgraded to HettInfo II. To capture a complete centrifugation process, we recommend the Add-on Kit plus Scan Module.

* Older models on request.

Articles	Includes
Add-on Kit for model 5005	Firmware, touch display, front panel
Add-on Kit for model 5005-80	Firmware, touch display, front panel
Add-on Kit for model 5005-90	Firmware, touch display, front panel

Cat. No.
E4378
E4393
E4416

Learn more: hettweb.com/hettinfo

Contact your local Hettich representative:

© 2019 Hettich Instruments, LP. All rights reserved. Lit. No. PS-HETTINO-EN-09-2019.

Hettich products are designed to help you achieve optimal results for your application and are built to perform to the specifications outlined in the operators manual. For application specific information and settings, please refer to your organizations standard operating procedure. As always, our Hettich representatives are here to help determine which Hettich products and/or accessories best fit your laboratory requirements.